


# gwrappon PF-1 User's Manual

## With Warranty

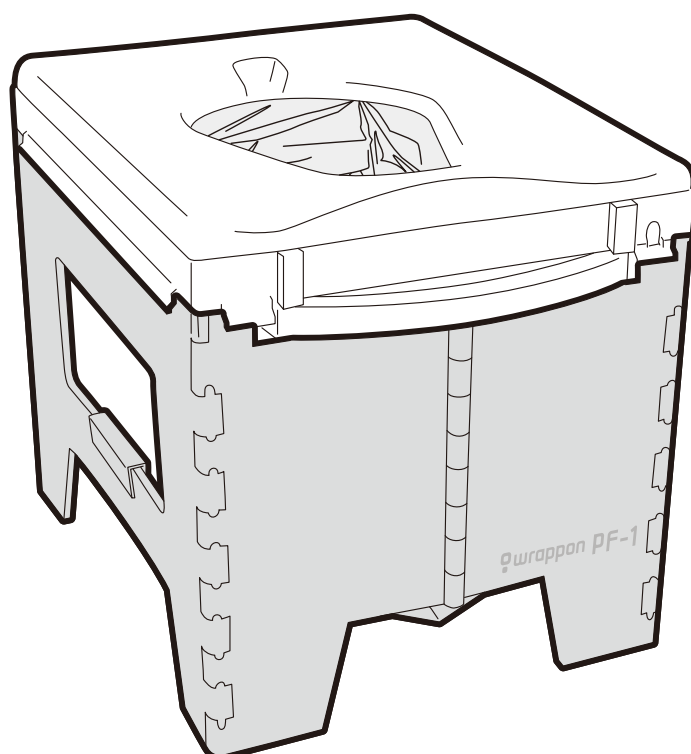

Thank you for purchasing Wrappon PF-1.  
Please read this manual carefully and use it correctly.  
After reading, please keep it in a safe place.

Product name	Wrappon PF-1
Product code	PF1SE001JH(Navy)
Rated voltage	AC100V
Rated frequency	50/60Hz
Rated power consumption	35W
Size	Width 376×Depth 510×Height 400mm
Product weight	approx.4.2kg
Material	Polypropylene
Maximum load capacity	100kg

## Table of Contents

### Before use

Safety precautions	.....2
Checking parts and accessory	..4

### How to use

How to assemble	.....5
How to set a garbage bag	...7
How to use	.....8
Handling care	.....9






### Maintenance and storage ·10











### Q&A .....11

### About warranty .....12







# Safety precautions

Before using the product, read this User's manual thoroughly and to ensure correct use. The precautions described here represent safety related issues and must be observed.




Precautions display and meaning		Contents to be protected is divided and explained by pictorial display.	
 <b>Warning</b>	Using the product incorrectly and ignoring this symbol may cause death or serious injury.		This type of symbol means "PROHIBITED". The one shown on the left means "Do not disassemble".
			This type of symbol means "REQUIRED.". Left sign means "Must be Observed".
 <b>Caution</b>	Using the product incorrectly by ignoring this symbol may cause injury or property damage.		This type of symbol means to enforce actions based on instructions. "Unplug the power plug" is shown.

 <b>Warning</b>		
 <b>Must be observed</b> <b>AC adapter must insert into AC100V outlet.</b> If the insert is incomplete, it causes a fire caused by electric shock or overheating.	 <b>Must be observed</b> <b>Do not use AC adapter other than specified one.</b> If used different one, it causes machine damage and is dangerous.	 <b>Must be observed</b> <b>Use it in a horizontal and stable place.</b> If you set it unstable place, you may fall and cause an injury.
 <b>Wet hands prohibited</b> <b>Do not touch AC adapter or switch with wet hands</b> This causes electric shock.	 <b>Wet prohibited</b> <b>Do not splash water to the product directly.</b> This causes an electric shock.	 <b>Prohibited</b> <b>When thunder strikes, do not touch the body or AC adapter.</b> This causes an electric shock.
 <b>Must be observed</b> <b>Do not put coagulant into your mouth.</b> Coagulant is not food and drink. For detail, please follow the precautions for using coagulant.	 <b>Disassembly prohibited</b> <b>Do not disassemble or repair and modify.</b> High voltage is applied inside the product so these cause electrical shock or fire and these may cause injury due to malfunction.	 <b>Unplug</b> <b>When not used for a long time, please unplug AC adapter from the outlet.</b> It may cause an electric shock or electric leakage fire due to insulation deterioration.


## ⚠ Warning

<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Do not damage codes and plugs.</b></p> <p>Do not scratch. Do not modify. Keep away from high temperature. Do not bend forcibly. Do not twist, nor pull. Do not put on the heavy item. Do not bundle.</p> <ul style="list-style-type: none"> <li>● If you use it while it is damaged, it may cause an electric shock or short circuit.</li> <li>● When unplugging, please hold the plug and pull it straight out.</li> <li>● For power cord, etc. repairs, please ask the store where you purchased the product.</li> </ul>	<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Do not overload the product.</b></p> <ul style="list-style-type: none"> <li>● Do not stand on the product.</li> <li>● Do not move with people on it.</li> <li>● Do not use with more than 2 people at the same time.</li> <li>● Do not put heavy objects on.</li> <li>● Do not give a strong impact.</li> </ul> <p>The product may be damaged and cause injury as you fell down.</p>	<p style="text-align: center;"> <b>Must be observed</b></p> <p><b>Clean the dust of AC adapter regularly.</b></p> <p>When dust collects, insulation becomes poor so it heats up and causes a fire. Unplug AC adapter and please clean the dust once in a month.</p>
<p style="text-align: center;"> <b>Fire-prohibited</b></p> <ul style="list-style-type: none"> <li>• <b>Do not install near fire</b></li> <li>• <b>Do not throw flammable</b></li> </ul> <ul style="list-style-type: none"> <li>● These cause short circuit or ignition and leakage of human waste due to product deformation so please be careful for at near the stove.</li> </ul>	<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Those who need assistance in moving and children, do not use it alone.</b></p> <ul style="list-style-type: none"> <li>● When you play a prank, it causes unexpected injuries.</li> <li>● When the product falls, it may cause injury.</li> </ul>	<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Do not use if you weigh more than 100kg.</b></p> <p>The product may be damaged and cause an injury.</p>


## ⚠ Caution

<p style="text-align: center;"> <b>Must be observed</b></p> <p><b>Before use, make sure to check for the presence of waste bag.</b></p> <p>When used without waste bag, leakage of human waste may cause product failure.</p>	<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Do not use anything other than a dedicated waste bag.</b></p> <p>If used other bags, it cannot be sealed or cause leakage and product failure.</p>	<p style="text-align: center;"> <b>Prohibited</b></p> <p><b>Do not leave it wet or dirty.</b></p> <p>Spots and mold occur so the product may damage.</p>
--	---	---


# Checking parts and accessories

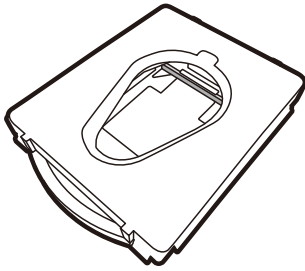


Toilet seat unit  
1 piece

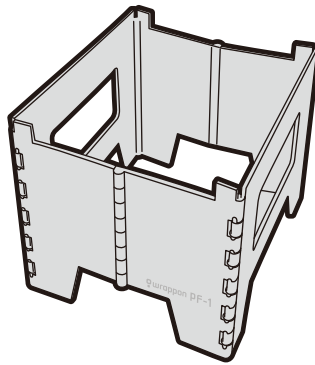


Main Body 1 piece  
Fold and store

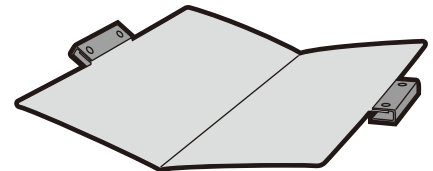

Pedestal 1 piece

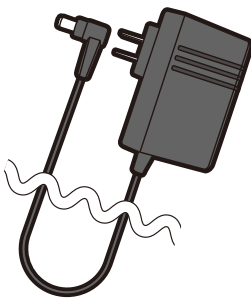


AC adapter 1 piece

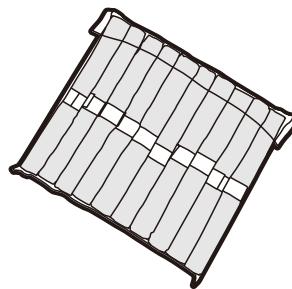

Dedicated waste bags  
3 bags(10 pieces per bag)



Dedicated coagulant  
30 pieces per bag

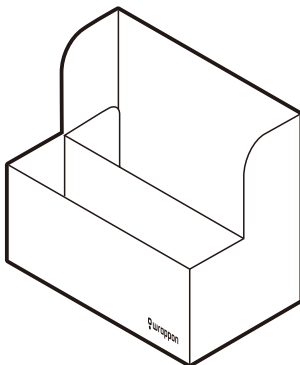

Consumable stand 1 piece


※Please use this as storage for coagulant and waste bags

# How to assemble




For more detail, please see the video of how to assemble of QR code on the left.

- 1 Spread the main body
- 2 Put the pedestal into the main body.

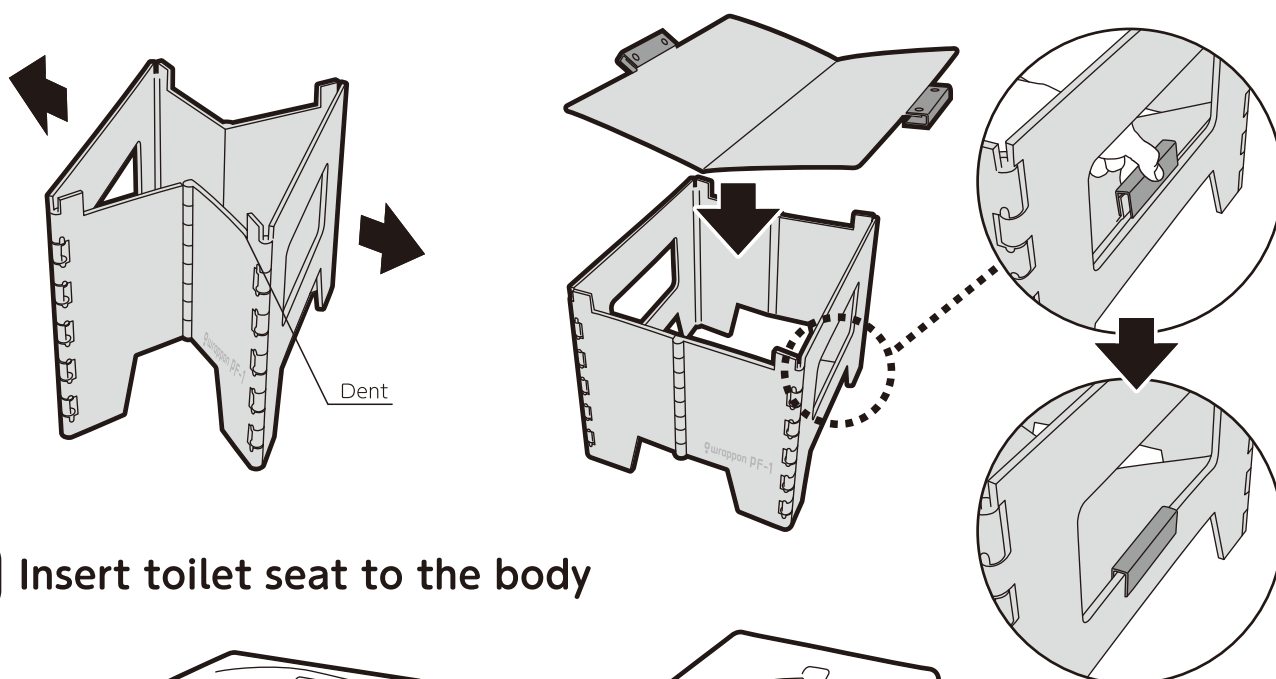

- 3 Insert toilet seat to the body

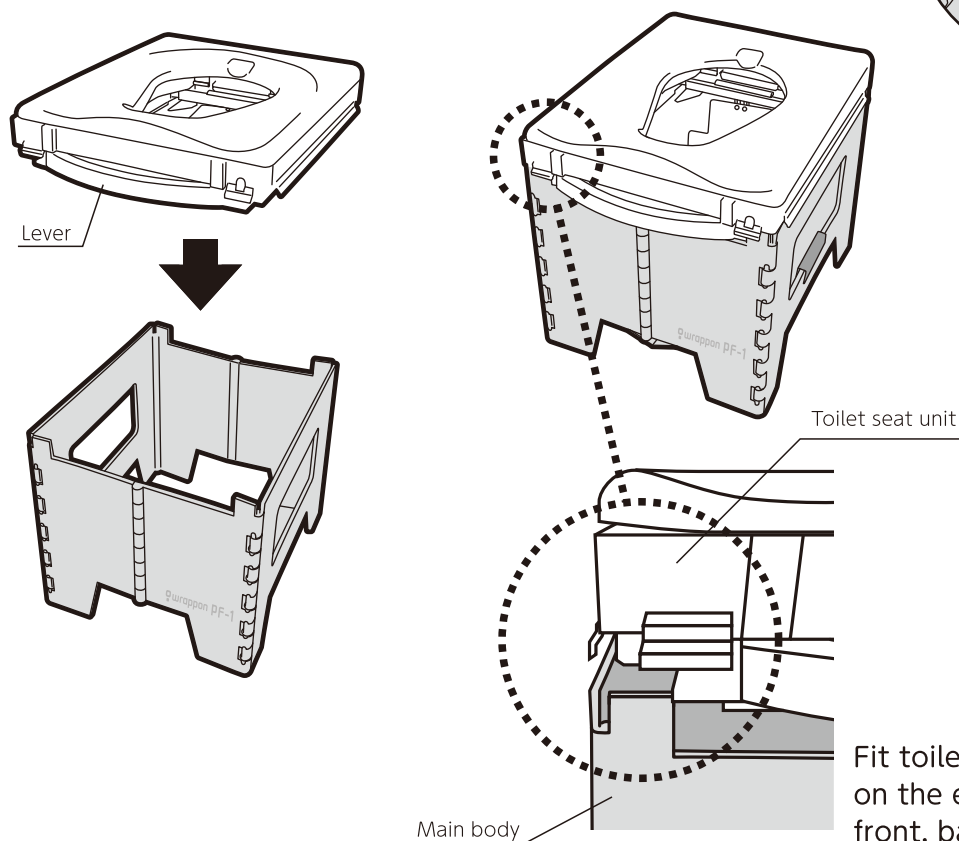


Fit toilet seat unit in firmly on the edge of the body front, back, left and right.

※When used due to insufficient fitting, it may cause a fall or injury.


# How to assemble(continued)

## 4 Connect AC adapter

① Please insert AC adapter plug firmly into the main body.

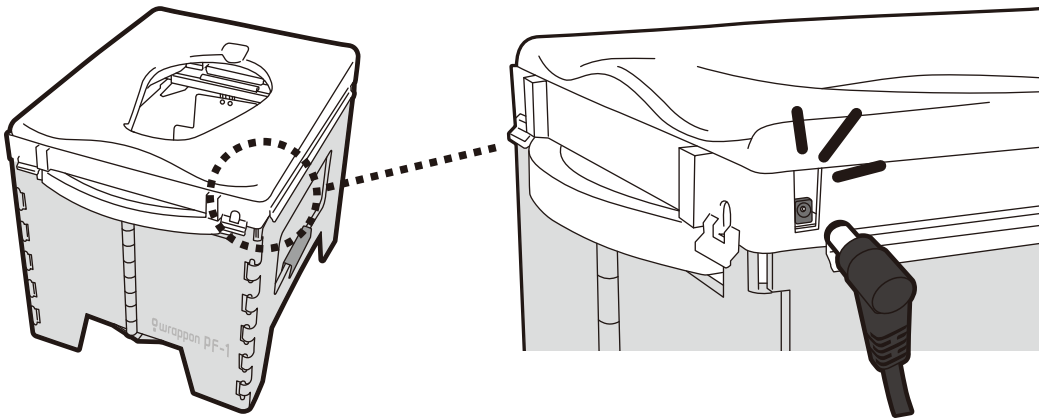

② Please insert AC adapter plug firmly into the outlet.

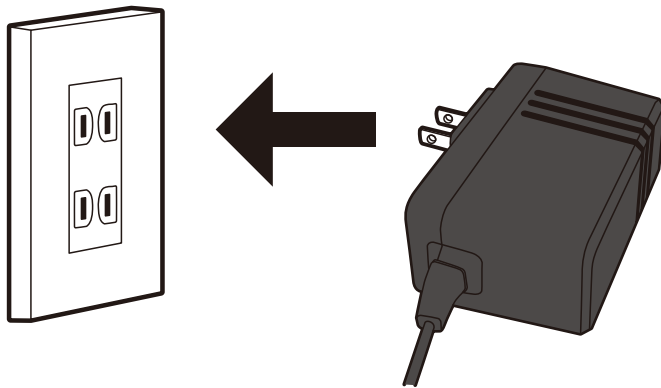

※ When using a battery, please insert the battery plug.

# How to set the waste bag





For more detail, please see the video of how to use and assemble on the attached QR code on the left.




## 1 Lift the toilet seat

Push the lever until the line displayed on the lever is hidden.

[View of the lever part from the above]

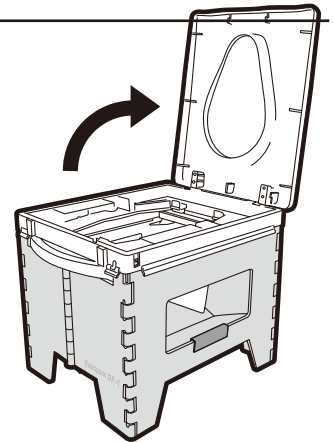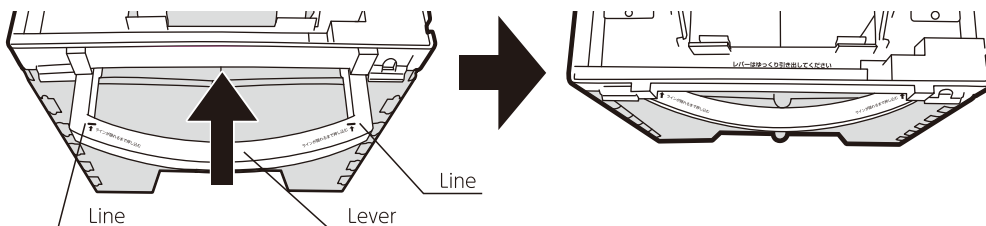

## 2 Open a waste bag into a hexagon

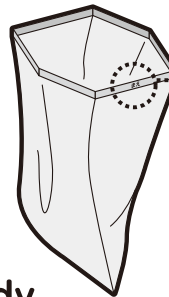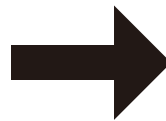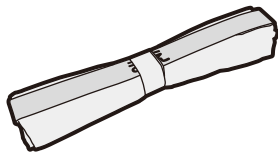

まえ  
Front

## 3 Set the waste bag on the main body

※When lowering the toilet seat, please be careful not to get your hands caught.

[View from the above]



[Clip part]

Back  
うしろ

Waste bag


Front  
まえ

The front of waste bag becomes the front of the main body(Lever side).

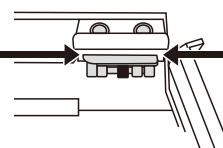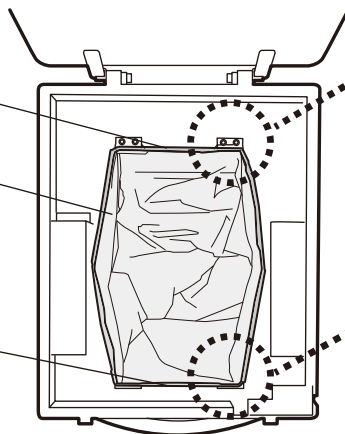

Insert the cardboard part of the waste bag into the clip.

[Groove part]

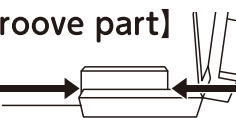


Insert the cardboard part of the waste bag into the groove part.

## Caution when setting

As shown in the figure on the left, when the left and right of the waste bag are installed inside the toilet, push back by hand to eliminate the gap. If there is a gap, human waste can leak from the waste bag.

Push back to eliminate the gap

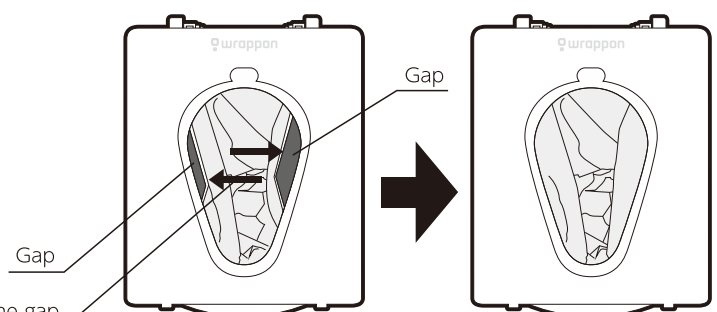



# How to use



For more detail, please see the video of how to use on the attached QR code on the left.

- When using, insert a plug into the outlet. ● Please sit in the posture of a Western-style.
- Do not stand and urinate. For urination, please sit down.


## 1 Pour the dedicated coagulant into the wast bag.

※You can also throw in the coagulant bag as well.



## 2 Relief yourself.


Sit on the correct position on the toilet seat.

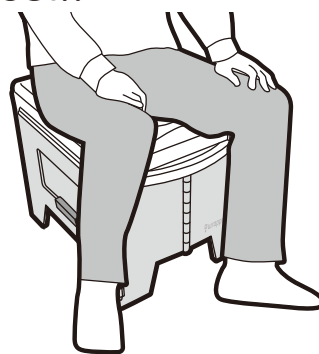


## 3 Lift up the toilet seat. While holding the main body by one hand, pull out the bar and up keep folding up.

※Pull out the lever slowly.(Pull out in about 3 seconds.)

When welding starts, you can hear beep sound and green light flashes. When completed, you can hear three beep sounds, then, light is lit to complete the process.

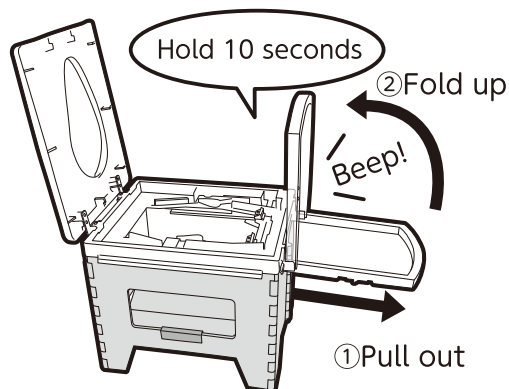


【Lever folded up state(View from the side)】

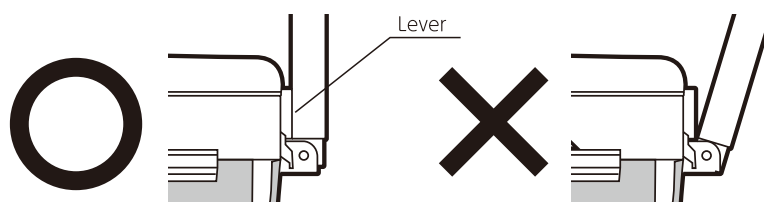


If you release the lever during welding, beep error sound is heard and light flashes fast. In that case, keep the lever holding up and keep welding.

## 4 Push the lever slowly and take the waste bag

Please dispose of the removed waste bag according to the sorting rules by each municipality.

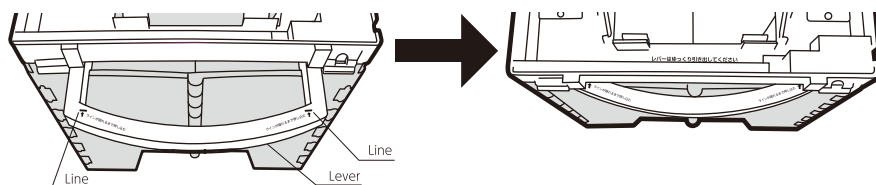

Push the lever until the line displayed on the lever is hidden.

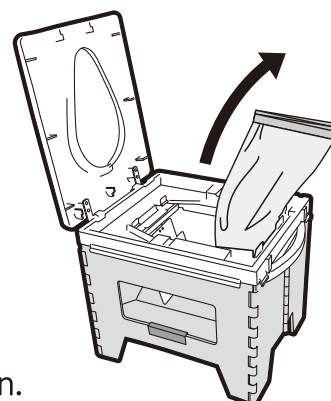




# Handle with care

## About the lever operation during compressed process.

Please pull out the lever slowly. (pull out in about 3 seconds)  
If you pull it too quickly, it may cause processing malfunction or more air in the waste bag.

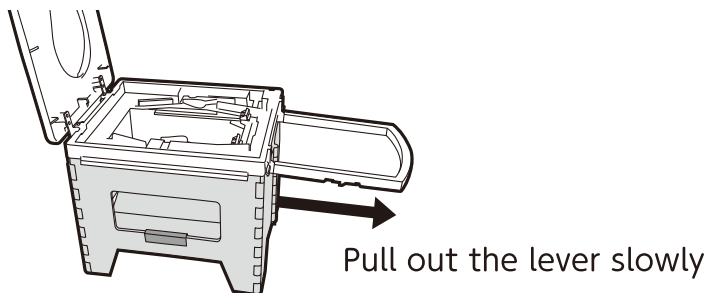


Pull out slowly,  
then, less air  
in the bag.




Pull out quickly,  
then, more air  
in the bag.

## About thermocompression process

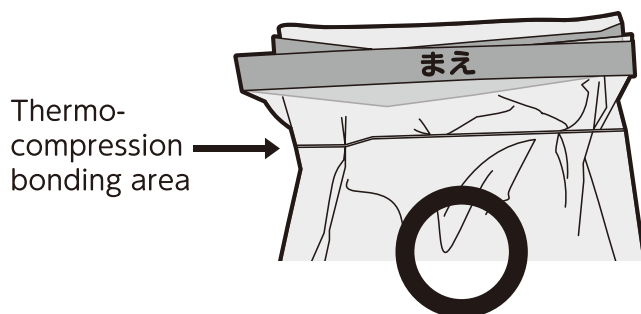


【Normal thermocompression】

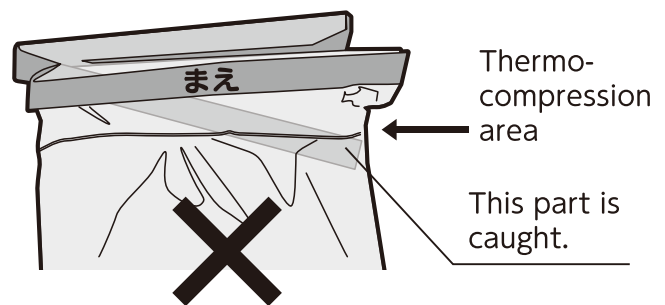


【Thermocompression area is caught in the part so it is incomplete.】

## Re-thermocompression

If thermocompression is incomplete (air leaks), please reprocess the removed bag as shown in the figure below.

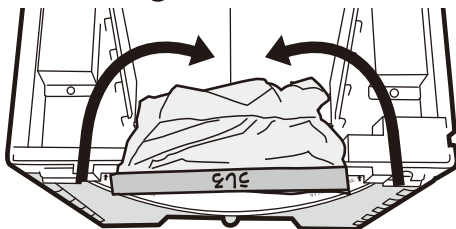

Re-thermocompression bag is hanging down as shown in the above figure and pull out the lever. Then, keep hold it as folding up.

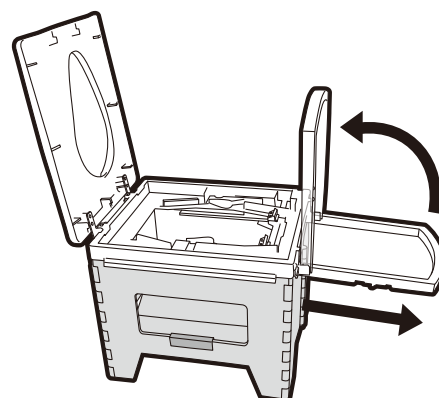

When thermocompression starts, beep is heard and green light flashes.  
When finished, three beeps sound is heard.  
Light is lit to complete the process.



For more detail, please see the video of how to use and assemble on the attached QR code on the left.

# Maintenance and storage

## About maintenance


- At the time of maintenance and inspection, make sure to unplug AC adapter from the outlet.
- Clean the dust of AC adapter regularly.
- Do not splash water directly on the main body. (The toilet seat can be washed with water.)
- When cleaning, please use water or a neutral detergent.
  - ※ Alcohol (disinfectant alcohol) and hypochlorous acid can be used.
  - Do not use toilet detergent, benzene, thinner, cleanser. etc.

## Cleaning a toilet seat


The toilet seat can be removed and washed.

### [How to remove the toilet seat]

- ① Remove the two hand screws behind the toilet seat.

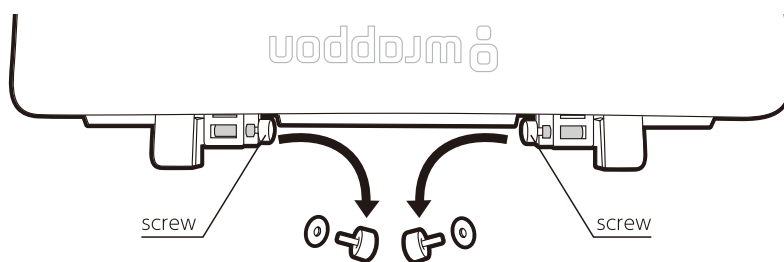

- ② Remove two shafts in the direction of the arrows, and remove the toilet seat.

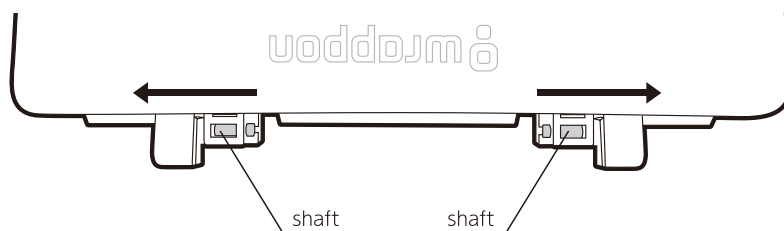

Installation is the reverse of removal.

※ After cleaning, make sure you wipe off any water on the toilet seat metal fittings. Remaining water can cause rust.

## About storage

Because it causes deterioration and damage due to falling so please do not use or keep in the following places.

- Places where it gets wet.
- Places where it is affected by air containing high temperature, humidity, direct sunlight and salt.
- Places where there are tilted, vibrated and get impact.

※ If you do not follow the correct usage such as the above precautions, we are not responsible for quality.

---

**Q1** What kind of storage location should I use?

**A1** Polypropylene materials used by the toilet seat unit, the main body and the pedestal deteriorate even if not used and even without using under the environment of hot temperature, humidity, direct sunlight and salt. Therefore, please avoid those places as a storage.

---

**Q2** What is the cause of poor thermocompression?  
If it becomes defective again, how to deal with?

**A2** Poor thermocompression happens when tissue paper or cardboard part of the waste bag are caught in the thermocompression area.  
When this happens, please refer re-thermocompression at page 9 to make re-thermocompression of the waste bag again.

---

**Q3** Can I make an additional purchase for the consumable products?  
Q3

**A3** Yes, you can.  
At the time of purchase of the product, 30 sets of the waste bags and dedicated coagulant come with the product. But, you need to purchase more as the number has decreased due to use and if there is a shortage for the estimated number of people.

---

**Q4** Can you repair the toilet seat unit?

**A4** If the seal mechanism is damaged, we will deal with either repair or replacement.  
Please refer the warranty at the page 11 if it is for free repair or with fee.

---

**Q5** If the main body or the pedestal is broken, can I buy only damaged parts?

**A5** Yes, you can buy each part.

---

**Q6** Has the waste bag a capacity limit?

**A6** There are individual differences for the quantity of human waste so it is recommended to make thermocompression for each excretion.

---

**Q7** What should I do when the power cannot be taken or when the battery runs out during thermocompression?

**A7** Remove the cardboard part of the waste bag and tie it by hand.

---

# About warranty

<b>Warranty</b>			
<b>Customer</b>		Purchase date	
Name		Warranty period	
		One year from the purchase date	
Address		<b>Repair record</b>	
Telephone		MM/DD/YYYY	Details of repair
			Manager

※If there is an omission, please note that we cannot guarantee.

This warranty promises to provide free send back repair which is within the range specified in the following free repair regulations and within the warranty period.

※Send back repair is a repair method in which you request repairs and send the product. The product will be sent after repair. By this warranty this does not limit customer's right. Please note that we may use the user registration details which are registered by the customer for the service activities during the warranty period and safety inspection activities

## **[Free repair regulations]**

1. During the warranty period, in normal use according to the precautionary statement by the user's manual, attached levels, etc. to this product, then, a failure occurs in this product. Then, how to respond? According to the regulations, we will repair with free of charge. (Based on our judgment, including replacement of all or part of this product. (The same shall apply hereafter.) Before requesting repair, please make sure to attach this warranty card to the product.

※Note that this warranty will not be reissued so please keep it in a safe place.

2. If you wish the product repair, after registering as a user, please contact the customer service help desk (user mail support) ※ This product support is by email only.

3. Even within the warranty period, in the following cases, these will be repaired with a fee.

1) If not registered.

2) No warranty presented

3) No serial number, no date of purchase, no customer's name, rewritten words

4) Failure or damage due to transportation, move, drop and other improper handling after purchase.

5) Failure or damage due to fire, earthquake, salt damage, gas damage, wind and flood damage, lightning strike, abnormal voltage and natural disaster

6) When the cause of failure is other than this product.

7) Failure or damage due to improper repair or modification by other than our company.

8) Failure, damage or dirt caused by improper handling, etc. which is different from the User's manual.

9) Although the cause of the failure is deterioration due to wear, exhaust, fatigue damage, direct sunlight and cleaning.

10) Failure or damage due to current (voltage, frequency) fluctuation other than specified.

4. This warranty applies Japan only.

Please be sure to check the contents. When using the product, it is considered that you have agreed to the contents described here.

## **About user registration**

Registration can be done from our website.

<https://wrappon.com/register>




## Before you call Wrappon safety call, please prepare the following information.

Please mention the following information of [Serial number] [Date of purchase] [Store you purchaed]

- ※ If it breaks down within the warranty period, we will conduct free repair if you show the warranty card. But, even within the warranty period, it is not subject to free repair regulations.
- ※ If you are not registered as a user, it is not covered by warranty.  
Before you call, please complete user registration.

### Serial Number attachment position

Please check the serial number which starts with N for 11 digit number.  
It is listed in the place on the inside left with the toilet seat etc.  
removed or the back side of the main body.

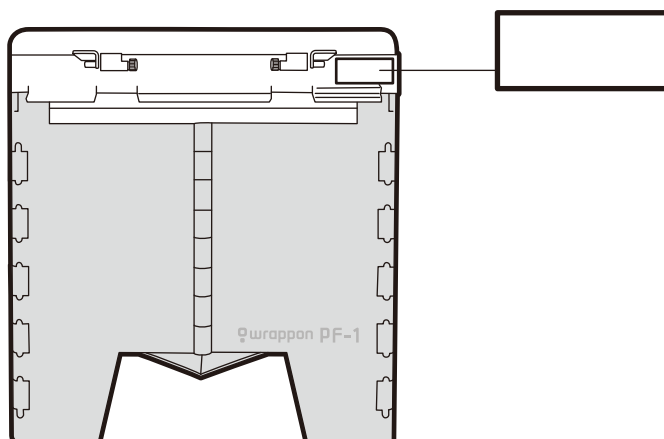

### Customer service center

About the product, we make absolutely sure to produce. But, if something goes wrong with you, please contact us after checking the contents thoroughly.

● User email support <https://wrappon.com/contact>





Ichibancho Tokyo Building  
21 Ichibancho Chiyoda-ku, Tokyo 102-0082 Japan  
Inquiry: [wrappon@nihonsafety.com](mailto:wrappon@nihonsafety.com)

※Please note that this use's manual may be modified for improvement without prior notice.

PF1SE001JH First edition in May.2022